

# Swing Out Centrifuge



**Brand:** HWLAB

**Product Code:** AC013HW

**Availability:** 2 - 3 Days

**Call for Price:** +2 01000 6546

## Short Description

HWLAB Swing-out centrifuge is widely used in physical and chemistry analysis, biochemistry...

## Description

HWLAB Swing-out centrifuge is widely used in physical and chemistry analysis, biochemistry, cellular and molecular biology, clinical Labs and blood donation centers. Horizontal separation offers many advantages over traditional fixed-angle separations, including tighter straight-line gel separations to reduce re-spins and re-mixes.


The unique swing-out rotor - streamline frontal area and high stability reduce the wind resistance effectively and load of motor, which is advantaged for reducing noise and temperature, make it quiet and protecting test samples.

## features

### Product features

- Maintenance-free brushless motor.
- LCD display of time and speed.

- Adjustable speed from 300r/min to 4000r/min.
- Timer, from 15s to 99min or continuous operation.
- For 24x 10ml/7m/5ml/3ml blood Collection tube.
- Lid automatically open at the end of the run to prevent sample warming and to reduce handing time.
- Three program button P1/P2/P3 for easy access to routine programs.
- The button rpm/rcf for easy to shift.
- Automatic rotor imbalance detection.

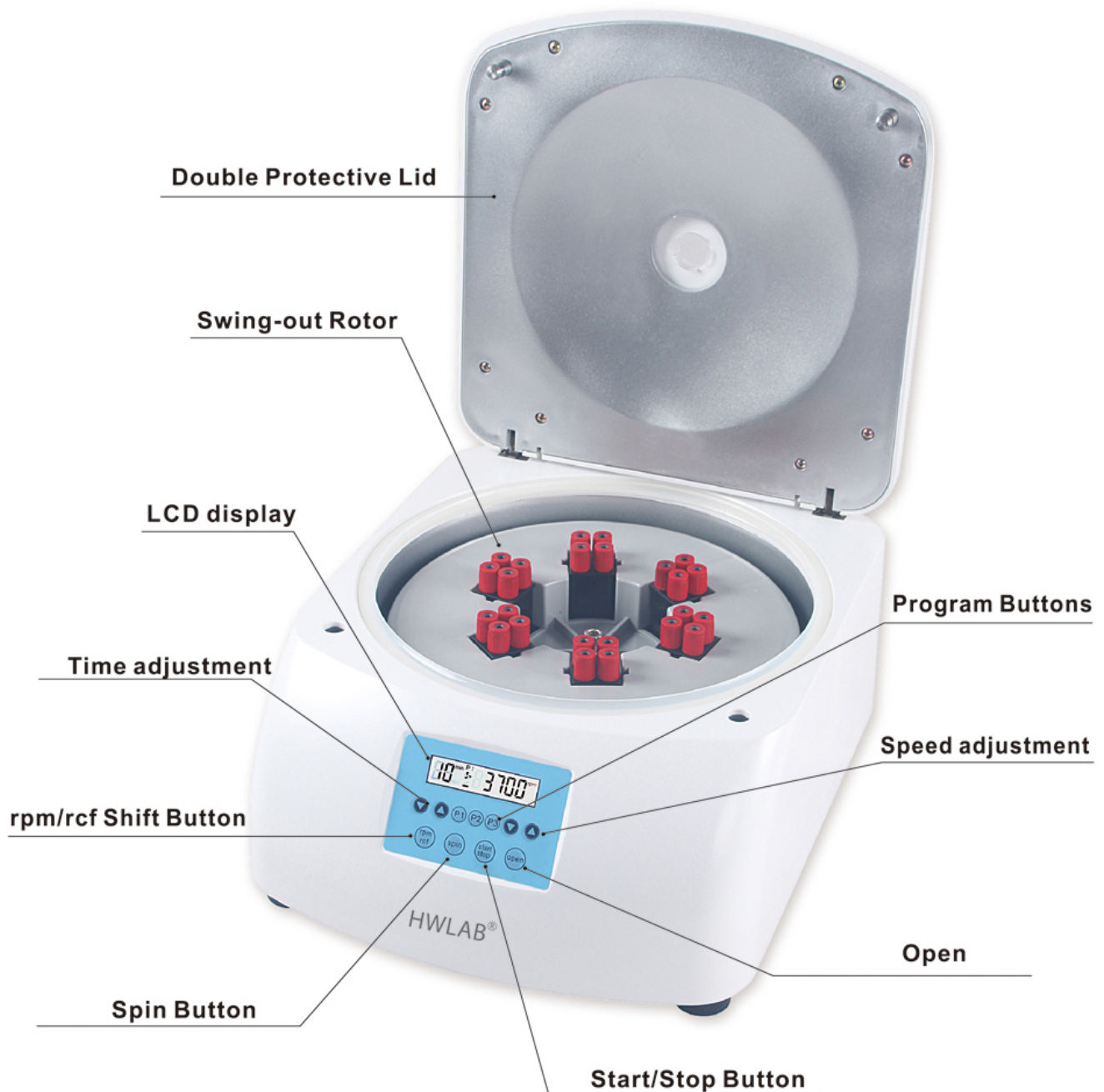

**Technical specifications**  
**Technical specifications**

**Maximum rcf**  
**Power supply**  
**Rotor capacity**

2600xg  
AC100V-240V/50-60Hz  
24x10ml/7ml/5ml/3ml blood collection tubes

**Speed**  
**Timer**  
**Motor type**  
**Dimension**  
**Weight**

From 300r/min to 4000r/min  
15s-99min / continuous operation  
Brushless motor  
450x375x225mm  
19KG

## Product Gallery

